



CS Voyager & AIO PC

Technical Training

October 13, 2021

© 2020 Carestream Dental LLC

Carestream Dental - Unrestricted Internal Use

Agenda

- **Technical Specification**
- **Installation**
- **Quick Start**
- **Maintenance**

Product Introduction

- ***CS Voyager -> Simply a mobile trolley***



- ***AIO PC -> All-in-one medical grade PC for intra-oral scanners***




CS Voyager + AIO PC



**Base diameter
575mm x 510mm**

**The holder of CS 3600/3700/3800
can be mounted on the handle**

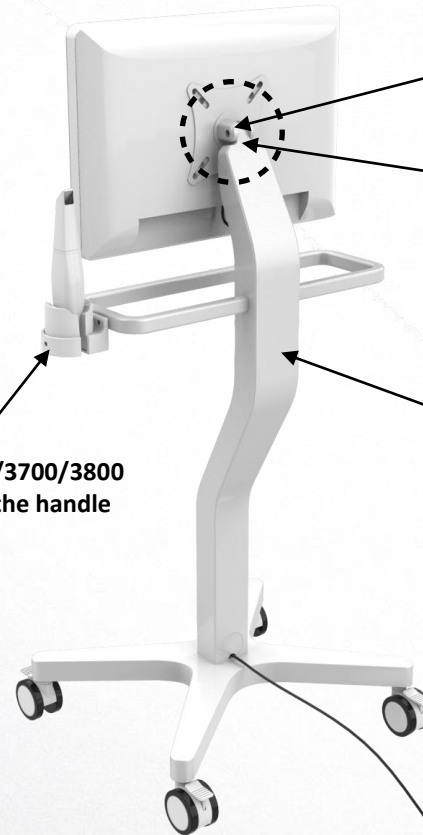

**AIO PC can be connected
by a standard VESA mount**

**Adjustable hinge
to adjust the
viewing angle
of the screen**

**Power Adapter of PC can be
hidden inside the stem**

**All the wheels can be locked
for stability**

Technical Specification

CS Voyager and AIO PC Dimensions

Item	Value
Dimension (L x W x H)	57.5 x 51 x 121 cm concerning the max. size. (without accessories)
VESA Mount dimension	Support 7.5x7.5 cm and 10x10 cm
Height(Center position of VESA mount to the ground)	116 cm
Height (Top surface of handle to the ground)	92 cm
Weight	13.7kg without accessories
Monitor up and down tilt range	Up and down to 50 degrees
Monitor side to side tilt range	Side to side to 50 degrees

!Note: The VESA mount is not able to be rotated.




AIO PC Technical Specification

AIO PC is 21.5" touch screen




Front View

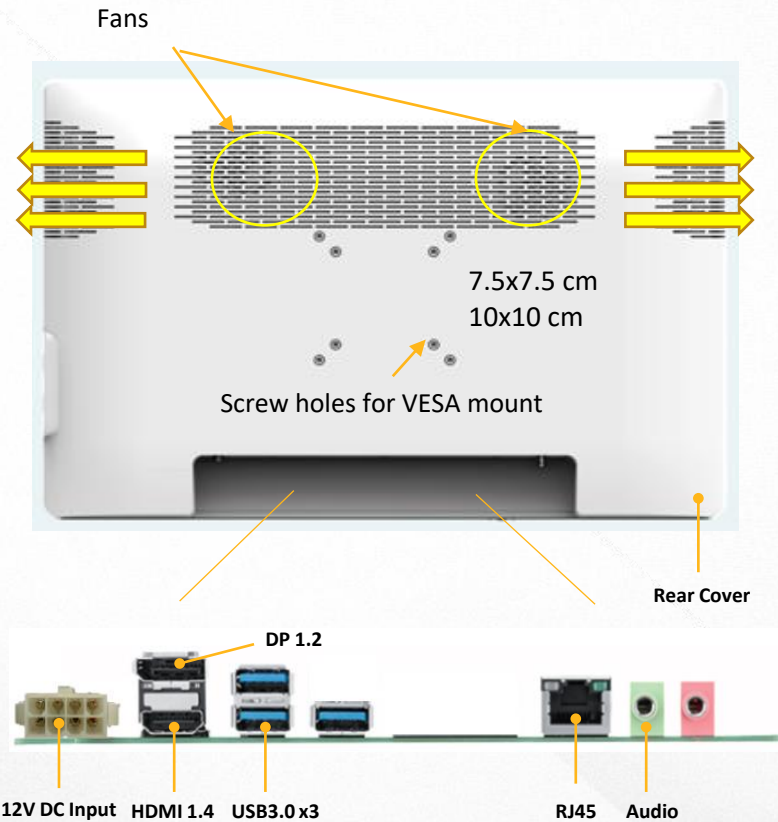

Rear View


AIO PC Technical Specification

Computing System	CPU	Intel Core TM i7-9700E Processor	
	Memory	16 GB DDR4-2666	
	Graphics	Intel UHD Graphics 630	
	Operating System	Win 10 IoT Ent LTSC	
	Storage	SSD 256GB (100 GB for C:\ and 156 GB for D:\)	
	WiFi+BT	IEEE 802.11 ac/a/b/g/n, Wi-Fi with Bluetooth 5.0 (M.2 Combo)	
	GPU	nVIDIA T1000 MXM GPU	
	Camera	Yes, 5M pixel camera	
	Microphone	Microphone 2V 400uA 114dB OD6X2.2 2P L200	
	Speaker	Speaker 2W 4Ω 78dB 32x23x14mm	
Display	Screen Size	21.5" TFT LCD (16:9)	
	Type	AHVN (IPS-like)	
	Resolution	1920 x 1080 (HxV)	
	Viewing Angel(degree)	178x178	
	Brightness(nits)	250 nits	
	Contrast Ratio	1000:1	
	Display Colors(bits)	16.7M (6bits + A-FRC)	
I/O Ports	USB Port	5x USB 3.0 (Rear: 3, Side:2)	
	HDMI	1 (HDMI 1.4)	
	DP 1.2	1	
	Ethernet (RJ45)	1	
	Audio	2 (Mic-in, Line-out)	
Physical Characteristics	Weight	6.6kg	
	Dimension	524 x 357 x 57 mm	
Environment	Temperature	Operating: 0°C~35°C	
		Storage: -20°C~60°C	
	Humidity	10~95% @40°C (non-condensing)	

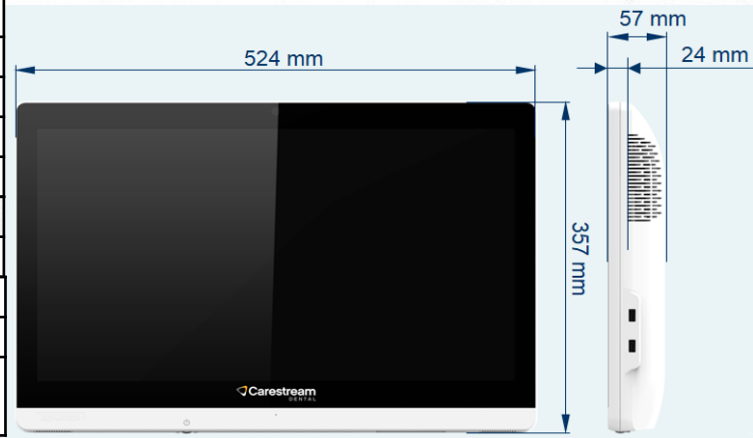


Installation

CS Voyager Components

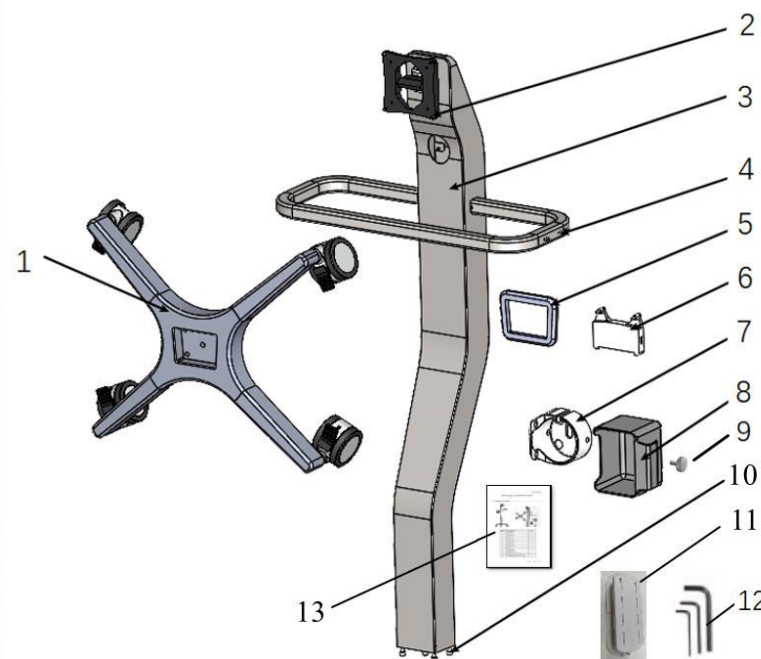

No.	Description	QTY.	No.	Description	QTY.	No.	Description	QTY.
1	Wheel base	1	6	Adapter base	1	11	Power strip hook	1
2	Screw M4x12	4	7	Charger holder	1	12	Installation tool (Allen wrench)	3
3	Display stem	1	8	Basket	1	13	Installation guide	1
4	Screw M4x20	4	9	Plastic knurled screw M4x10	4			
5	Silicone rubber ring	1	10	Screw M6x20	4			

AIO PC Standard Package

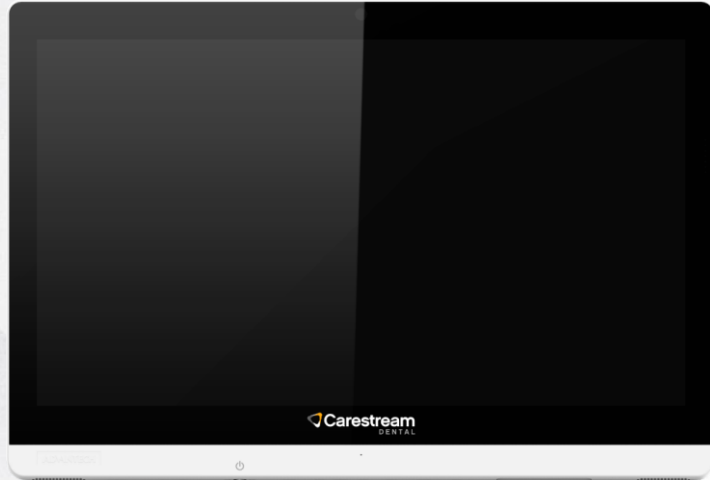

All in one PC



Power Adapter



USB Extension Cable (0.2m)
USB to Type-C Adapter

Input	100-240 V, 2.9-1.25A, 47-63 Hz
Output	220W/Max
Cable length	80 cm


***Power Cord is not packed with the PC and needs to be ordered separately**

CS Voyager Installation


- Display stem installation

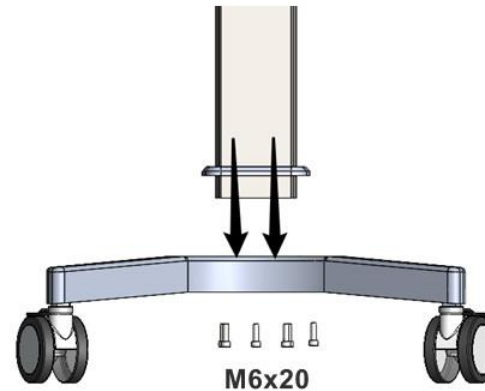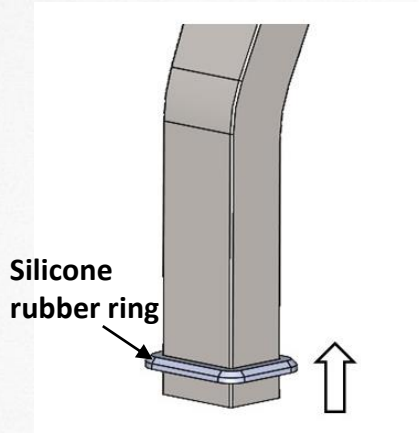

- Lay the stem down so you can install it easily and put foam beneath the cart to avoid scratches.

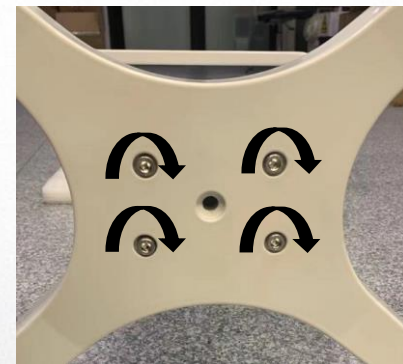

AIO PC Installation – With Trolley

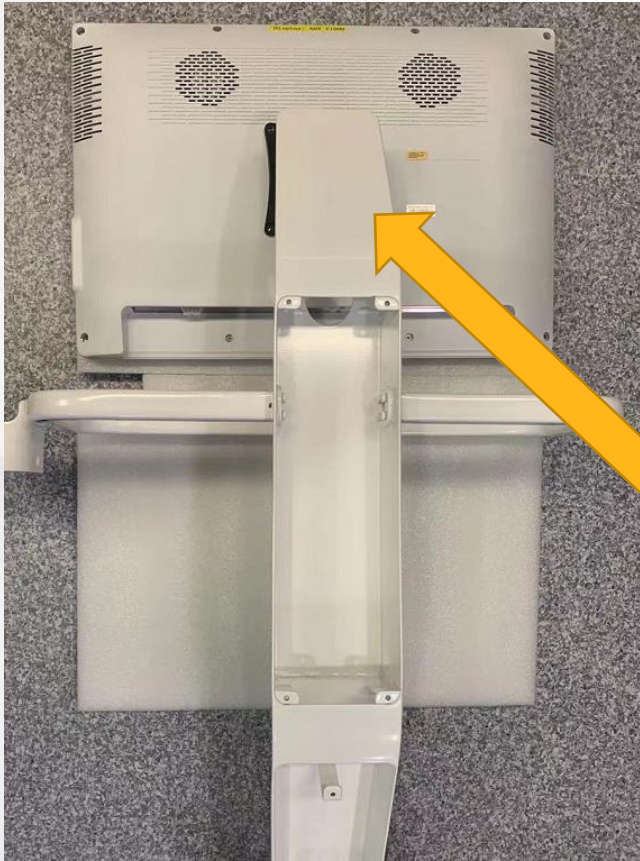

- Place foam on the ground and place the AIO PC on top of it.
- Gently laydown the trolley and align the VESA mount with the holes on the PC.
- Install and tighten 4 screws (M4 x 12) on the VESA Mount.
- Standup the trolley.





M4 x12


Accessory installation

➤ Adapter base

- ❖ For CS 3800 Wireless, mount the WiFi adapter. A USB 3.0 extension cord is required for the connection.
- ❖ For CS 3600/CS 3700, slide the power box into the adapter base

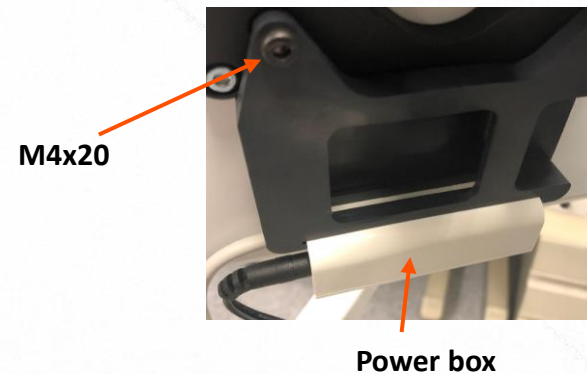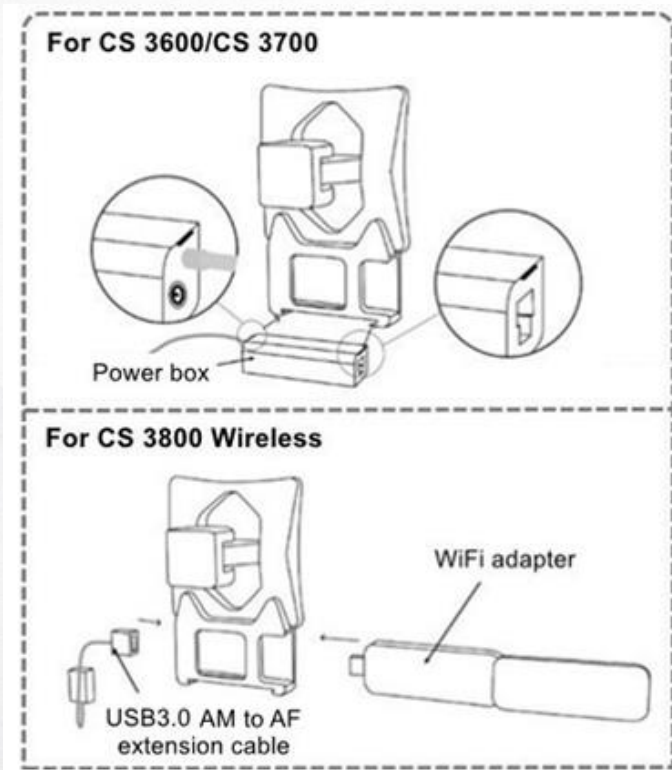


Accessory installation

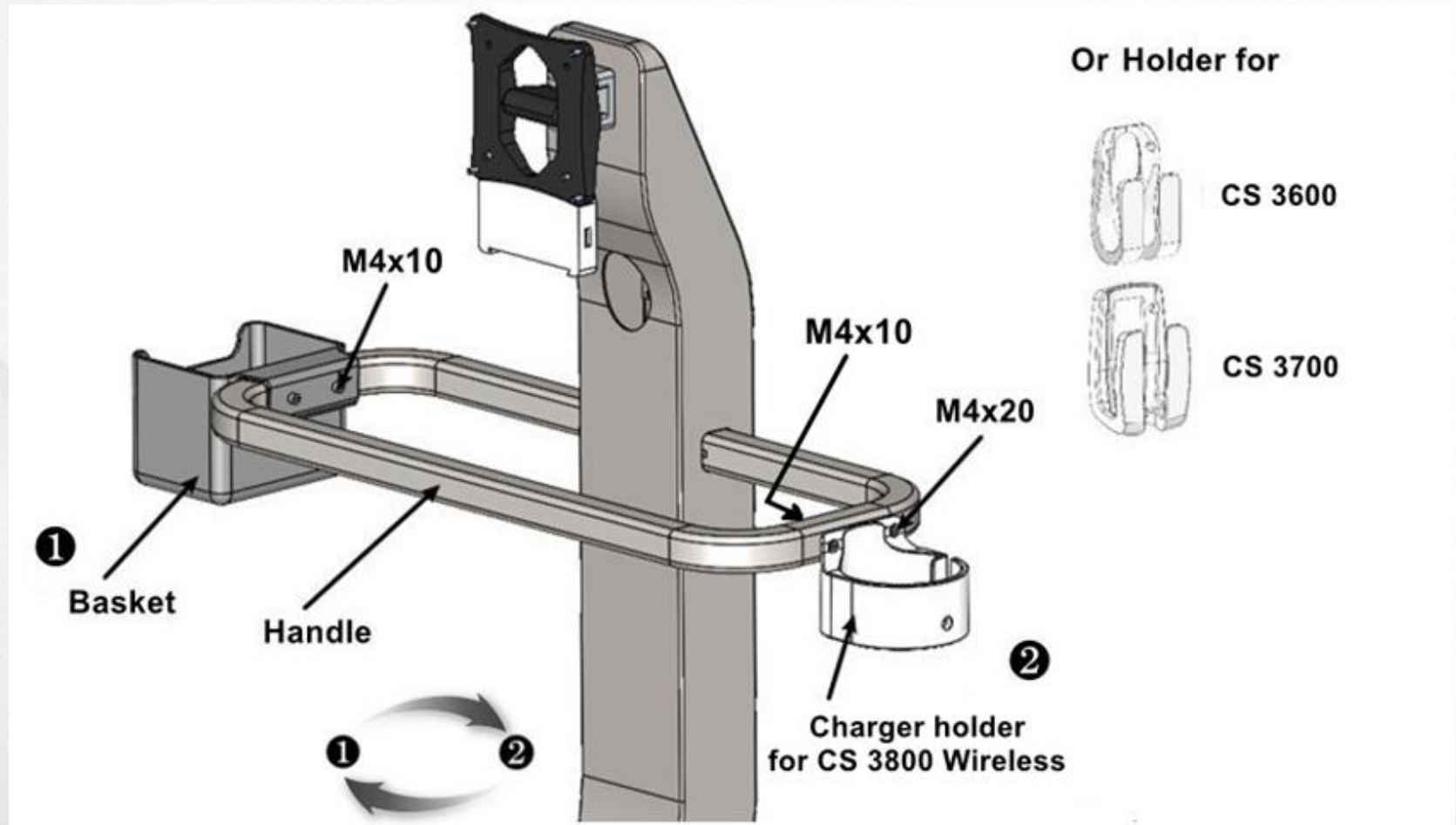


Cable storage and management

- AIO PC– Power adapter
- IO Scanner
 - CS 3800 Wireless – Power adapter for charge holder
 - CS 3600/CS 3700 – Power box and power adapter



AIO PC

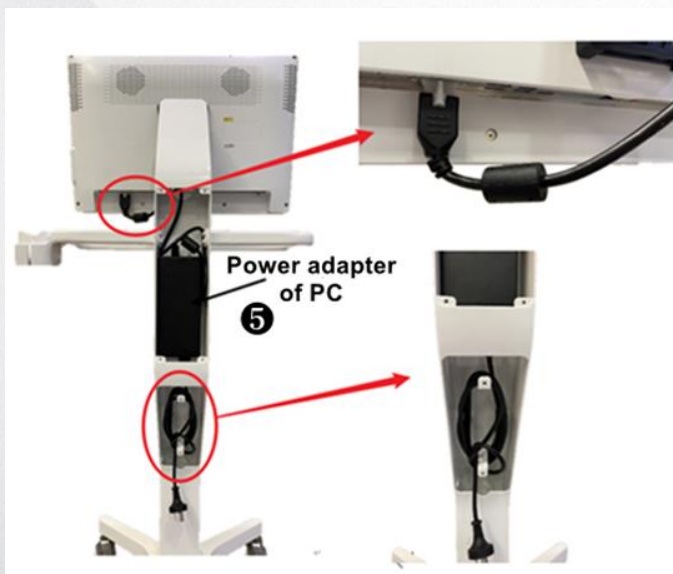

For CS 3800 Wireless

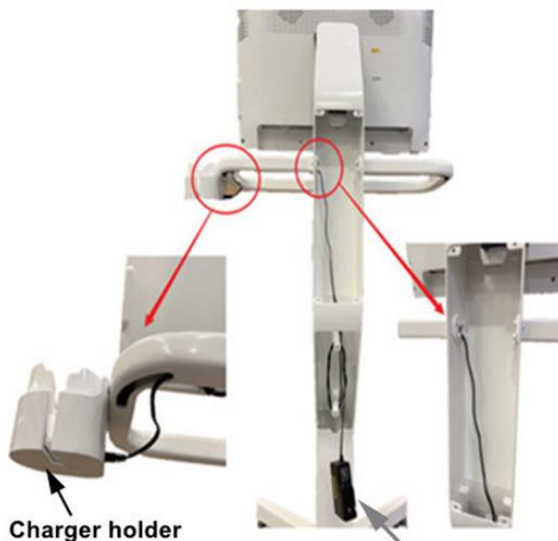


For CS 3600/CS 3700

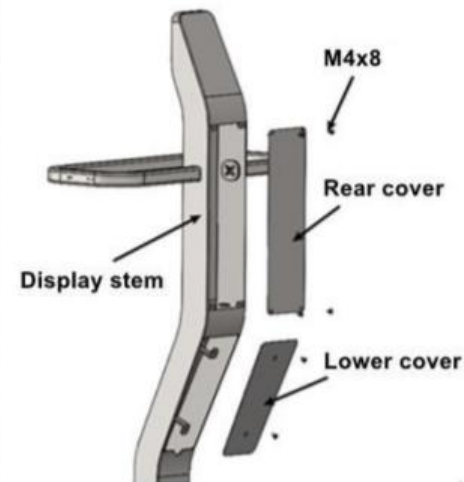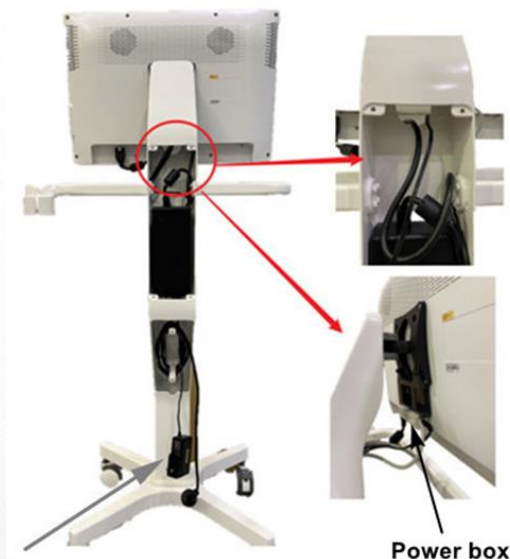


Cable storage and management

- Power strip hook installation
- Power strip cable storage

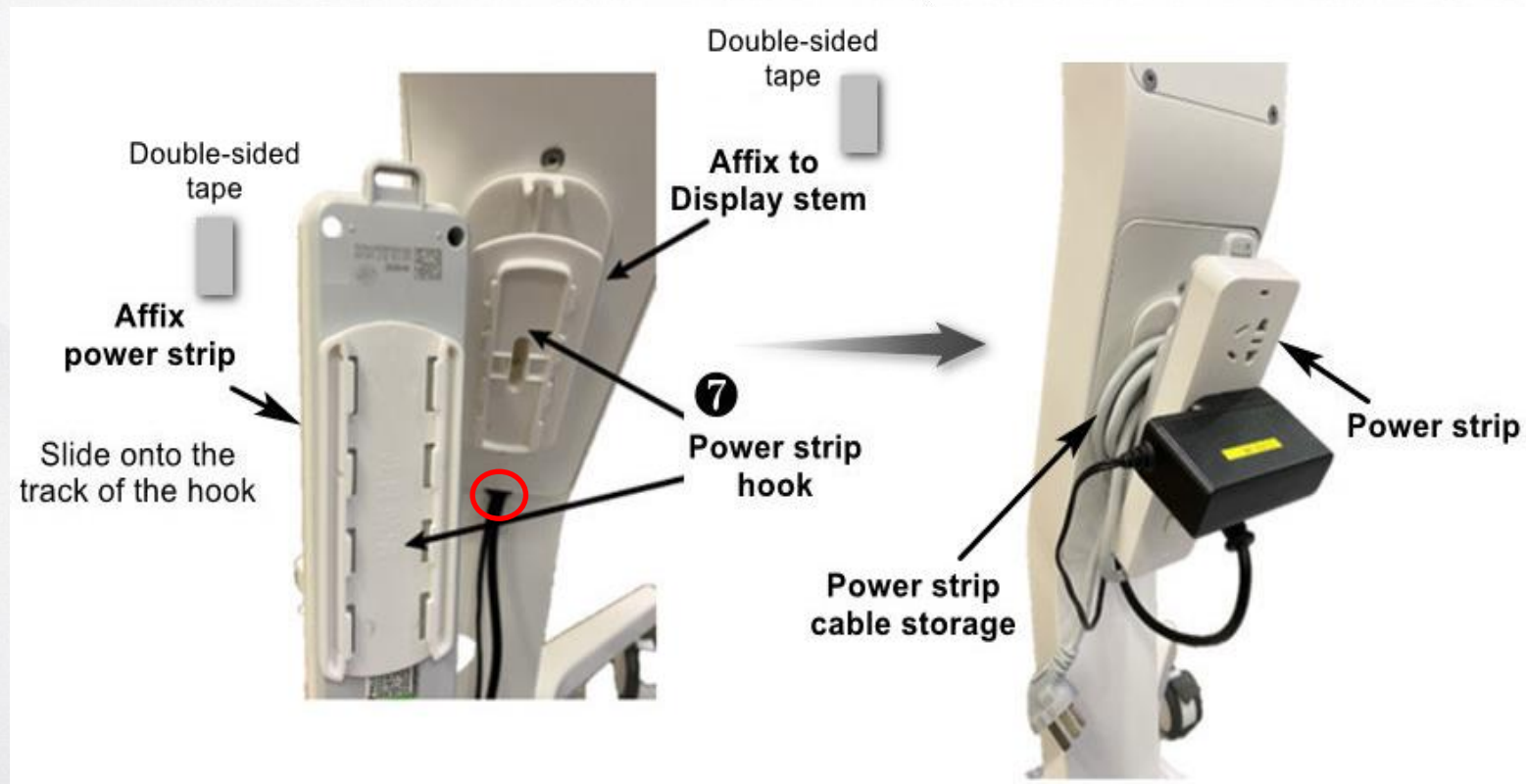


Two pieces of double-sided tape are included with the power strip hook.

Desktop Stand CAT# 5943782

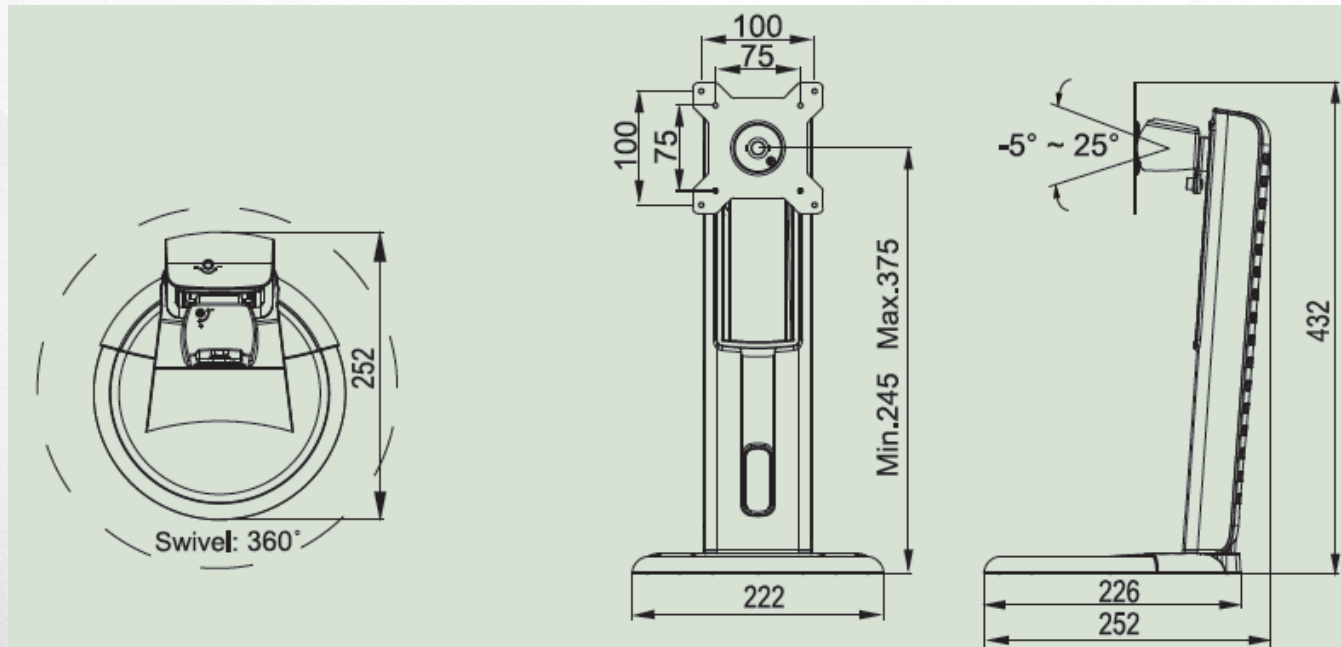


AIO PC Installation- With desktop stand

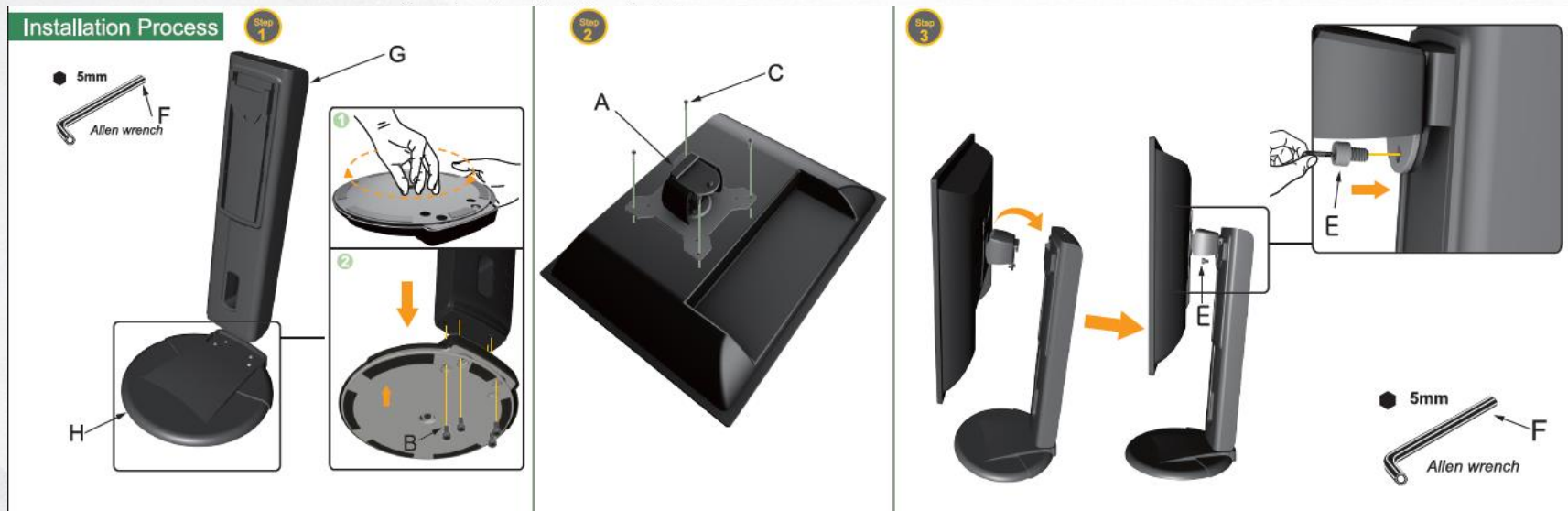


Swivel&Tilt / Rotate

Vertical Adjust

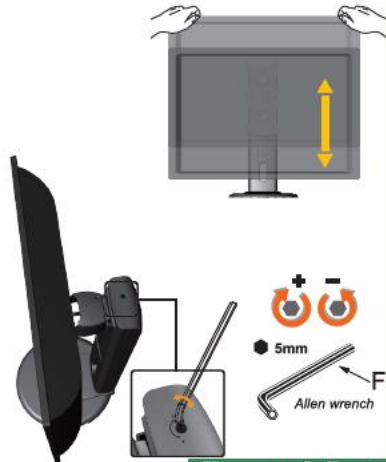


Torque Adjust

Tilt

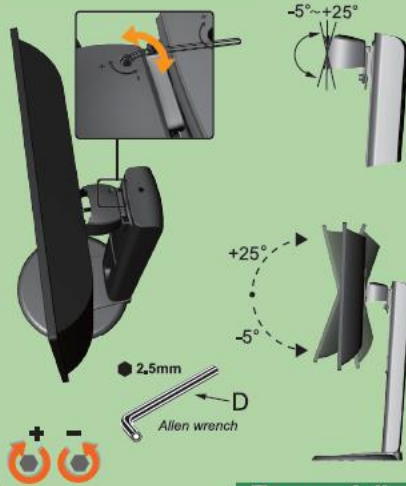

Torque Adjust

Rotate

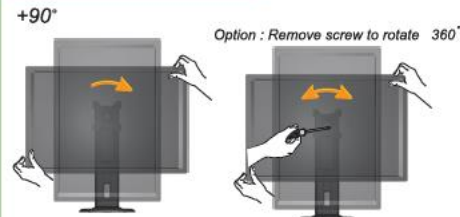

Swivel



Cable Management



Quick Start

AIO PC PC Setup

- Tear off the protective film on the screen before use to enable good touchscreen experience.
- Operation System: Windows 10 IoT Ent 2019 LTSC
- User account creation
- Software installation
 - CS Imaging software
 - CS ScanFlow
 - Any other software that is needed
- Database creation
 - CSI database
 - The folder to save exported image file and CSZX
- Backup
 - Please follow HT#126 to set up the backup of images

Recommend to set the database to another PC or sever.

Maintenance

Clean and Disinfection – AIO PC

Below wipes& solutions have been tested with the AIO PC.

1. AHP Accel TB
2. 75% Alcohol
3. Glutaraldehyde [C₅H₈O₂], 2% concentration
4. Isopropyl alcohol [(CH₃)₂CHOH], 70% concentration
5. Sodium hypochlorite [NaOCl], 10% concentration
6. CHLOR-CLEAN-1.7g NaDCC(1000 ppm available chlorine)
7. Clorox Healthcare® Hydrogen Peroxide
8. Clorox Healthcare® VersaSure
9. Clorox Healthcare® Bleach Germicidal Wipes

Clean and Disinfection- Trolley

1. 75% Alcohol
2. CaviCide
3. CaviWipes

Above solutions and wipes were tested by us.

!Note: Using untested solution may cause discoloration to the painting. Do not use the Sodium Hypochlorite which will cause yellowing.



Service Parts

- AIO PC Power Adapter, (SPCX761)
- Wheels of the trolley (SPCX762)
- Hinge assembly, SPCX763. Use of a 5mm Allen wrench is required



- Holder Kit, SPCX764

Includes the holder for CS 3800 and the adapter base

CAT# List

Cat #	Description
5943758	POWER CORD, EMEA
5943766	POWER CORD, US&C
5943774	POWER CORD, UK
6890313	POWER CORD, CHINA
6890321	POWER CORD, JAPAN
5943782	Stand for AIO PC
5943790	AIO PC All-in-One PC
5943832	CS VOYAGER
Consumables	
5944434	Basket, CS VOYAGER
Spare Parts	
SPCX761	AIO PC Power Adapter
SPCX762	Wheels Kit (4 wheels)
SPCX763	Hinge Assembly
SPCX764	Holder Kit



Thank You

